

natus<sup>®</sup>  
pediatrics

# CookieCap™

## CPAP Stabilizing Cap

Keeps CPAP Tubes Comfortably & Securely in Place

- Latex-free
- Can be placed over a hat
- Safe - no pins or rubber bands needed
- Easily adjustable
- Available in three sizes:
  - Micro-Premie**  
(18-25 cm head circumference)
  - Premie**  
(21-30 cm head circumference)
  - Full-Term**  
(26-34 cm head circumference)




COOKIECAP'S SOFT FOAM FABRIC HELPS ENSURE BABIES' COMFORT, AND KEEPS CONTINUOUS POSITIVE AIRWAY PRESSURE (CPAP) TUBES STABLE — WITHOUT HARMING FRAGILE NARES

### Keep CPAP dual tubes stable & ensure babies' comfort — with the unique CookieCap

- Gentle, secure attachment at temple reduces pressure on nares
- Chin strap prevents cap from sliding up and off head
- Chin strap stretches to allow nipple feeding & pacifiers
- Soft foam material is soft & stays in place
- Open top allows for scalp vein access & better bilirubin therapy
- Over-the-head strap is included as an option for additional stabilization

Your babies will love  
*CookieCap*—  
invented by the father  
of twin preemies

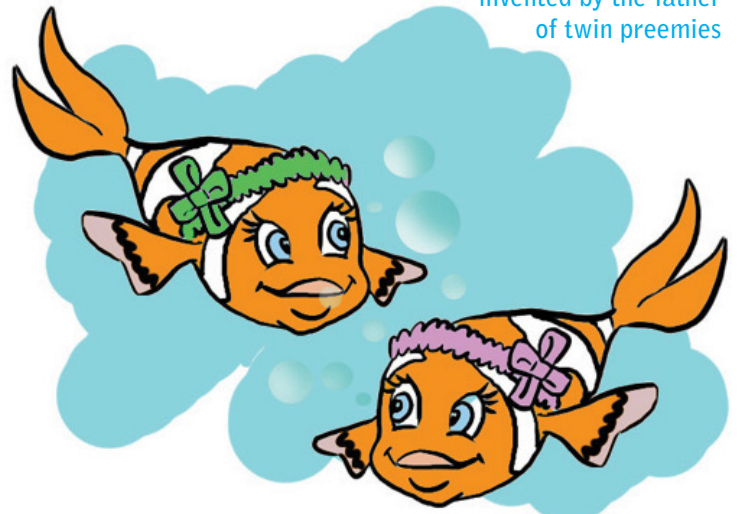

Hawaii Medical  
CookieCap™

Wish there was a safe, comfortable way to stabilize CPAP hoses and nasal prongs? There is – the Hawaii Medical CookieCap.



**TWO PIECES INCLUDED:**  
 • HEADBAND  
 • CHINSTRAP

FOAM HEADBAND WITH VELCRO® BRAND HOOK & LOOP – SECURES ATTACHMENT AT TEMPLE – REDUCES PRESSURE ON NARES



SOFT, COMFORTABLE CHIN STRAP PREVENTS CAP FROM SLIDING UP AND OFF HEAD, ALLOWS NIPPLE FEEDING AND PACIFIERS

CPAP CORRUGATED TUBES ARE GENTLY, YET SECURELY HELD IN PLACE VELCRO® BRAND HOOK & LOOP AND STRETCHABLE ELASTIC OVER-WRAP – NO SAFETY PINS OR RUBBER BANDS REQUIRED

PRODUCT SPECIFICATIONS:  
 Polyester polyurethane foam laminated with nylon looped fabric.  
 Latex-Free. Made in Mexico.

## Ordering Information



ITEM#	DESCRIPTION	QUANTITY
<b>For Dual Tube CPAP</b>		
1060021-15	Micro-Preemie Size (18-25 cm head circumference)	Box/15
1060022-15	Preemie Size (21-30 cm head circumference)	Box/15
1060023-15	Full-Term Size (26-34 cm head circumference)	Box/15
<b>Replacement Chin Straps</b>		
1060026-50	One size fits all	Box/50



The "Twin Cookies" Delaney & Madelyne 1 1/2 years old; old method of securing dual tubes with rubber bands and safety pins.



**Designed by Aaron Cook - father of twin preemies**

During his daughters' three-month NICU stay, Aaron watched Delaney and Madelyne struggle with CPAP. He designed the CookieCap to help stabilize CPAP tubes and reduce pressure on the nares. CookieCap helps maintain positive air pressure, while reducing septal breakdown and re-intubation.

Natus Pediatrics is proud to introduce CookieCap so all babies can comfortably benefit from CPAP. The innovative design helps ensure babies' comfort and keeps CPAP tubes stable – without harming fragile nares.

## We Welcome Your Feedback

As part of our ongoing focus on continuous improvement, we encourage you to share your thoughts on using our quality products. Please share your comments at [www.natus.com](http://www.natus.com) (click on Contact Us to send us your message).

### ALSO AVAILABLE:

#### Newborn Wrap Hats

The gentle fabric can be used to wrap around the baby's head to help prevent heat loss & still provide access to the scalp if required. The CookieCap stabilizing headband can be applied around the wrap hat to secure the dual tubes comfortably in place.



ITEM# OBM00043  
 5 x pink and 5 x blue, 10/pk



**natus**  
 pediatrics



Natus Medical Incorporated

1501 Industrial Road • San Carlos, CA 94070 USA | 1.800.303.0306 (US/CA) • Worldwide: 1.650.802.0400  
[www.natus.com](http://www.natus.com)