

natus[®]
pediatrics



LilyPad™

Pressure Relieving Memory Foam Mattress

- Replaces costly single use SheepSkin
- Economical:
Clean and re-use
Doesn't stain
- LilyPad meets small baby protocols:
Provides pressure relieving support
- Radio translucent overlay does not show on X-rays
- Available in **many** sizes and we'll custom make to your specifications



MEMORY FOAM DISTRIBUTES WEIGHT EVENLY AND REDUCES CAPILLARY BED COMPRESSION.
SOFT COVER IS WATERPROOF AND DOESN'T STAIN. LILYPAD COMPLIES WITH FEDERAL FLAMMABILITY STANDARDS.

The supple support and soft surface of *LilyPad™* keeps babies more comfortable and helps minimize sleep disruption. Cleanable cover doesn't stain!

- Pressure relieving mattress helps babies be more comfortable, reduces fidgeting and sleep disruption.
- Soft, responsive Memory Foam reduces pressure-induced soft tissue damage and prevents pressure sores in immobilized patients.
- Cleanable cover is treated with antimicrobial agent; Ultrafresh DM50. Cover passes anti-bacterial and anti-fungal standards.
- Economical – **can be cleaned and reused** between patients.
- Ideal for ECMO babies and high frequency ventilation.
- Latex-free.




A standard incubator mattress is like sleeping on a tabletop.
What do **YOUR** babies sleep on?

Hawaii Medical
LilyPad™



Occipital Region of Infant Mannequin Comparison of Interface Pressures

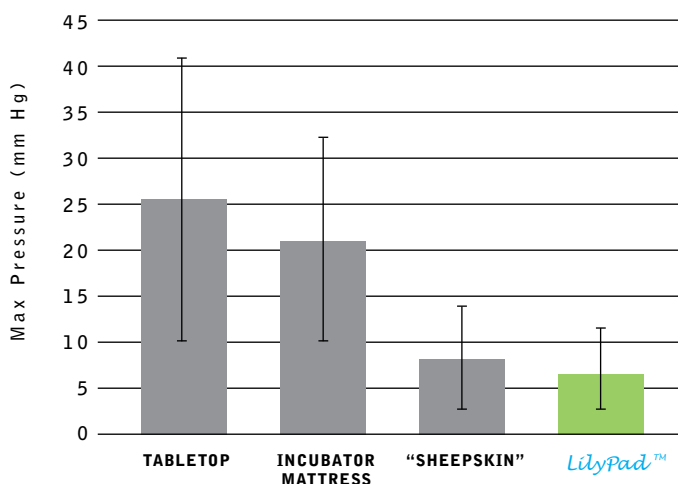

Support Surface	Interface Pressure (mmHG)	P VALUE
Incubator Mattress N=6	22+/-10.4	
TableTop N=6	26+/-15.2	0.41
"SheepSkin" N=6	8.6+/-5.5	0.0005
LilyPad™ N=7	7.1+/-4.5	0.0004
Feather Pillow N=8	4.6+/-3.3	0.0004

TESTS SHOW THE LILYPAD™ IS SIGNIFICANTLY SOFTER THEN A STANDARD INCUBATOR MATTRESS. THE REUSABLE LILYPAD™ OFFERS AN ECONOMICAL ALTERNATIVE TO SYNTHETIC "SHEEPSKIN" WITHOUT SACRIFICING PRESSURE RELIEVING PROPERTIES.

Now you can get a comfortable, pressure-relieving mattress designed for infants

Standard incubator mattresses are not significantly softer than a tabletop. Imagine trying to sleep on that for months at a time! It is no wonder some babies fidget. The low interface pressure of the LilyPad provides optimum pressure relief and maximum comfort for infants.

The fluid proof, anti-bacterial cover means you can clean and reuse the LilyPad between patients...saving you money while giving your babies the benefit of pressure-relieving support.

PRODUCT SPECIFICATIONS:
 Mattress – opencell polyurethane foam (Latex-free).
 Cover – polyester knit with urethane coating.
 DM50 as an antimicrobial agent. Antibacterial and antifungal.
 Cover Cleaning Instructions: Neutral soap and lukewarm water.
 See cleaning recommendations on web site.
 Cover Disinfection Instructions: Prepare and use dilute disinfectant and/or germicides in accordance with manufacturer's recommendations.

FLAMMABILITY TESTING:
 • Tested in accordance with ASTM E 1590-99.
 • Complies with the requirements of California Technical Bulletin 129.
 • Meets requirements of 16 CFR 1633 and CATB 603.

Ordering Information

Item #	Sizes in Inches 1.5 INCH THICKNESS	Equipment Application
WARMER MATTRESSES		
1030091	18.5 x 26	Air Shields® Open Bed Warmers Models 78-1; 78-2
1030092	18 x 29	Draeger® Babytherm® Open Bed Warmer Model 8004
1030093	21 x 26	Air Shields Resuscitaire® Radiant Warmer BRW & RW
1030095FP213	24.4 x 28.3	Fisher & Paykel® Open Bed Warmers Model 9001W213
1030095FP201	24.4 x 24.4	Fisher & Paykel Open Bed Warmers Model 9001W201
1030091A	18 x 24.5	Ohmeda® Open Bed Warmers (& Ohio®) Models 2001, 3300, 3400
1030096	23.25 x 29.5	Ohmeda Open Bed Warmers (& Ohio) Models 4300 and 4400
10300913	26 x 21 x 1/4	Acrylic Insert, standard size for warmers with xray opening. Sold individually.
INCUBATOR MATTRESSES		
1030094	13 x 26	Care Plus® Series and Care Plus® Access 3000/4000 Series
1030097	16 x 31	Air Shields Isolette C2000 & C2000E
1030099	18.75 x 24.75	Ohmeda Giraffe OmniBed (with rounded corners)
1030091A	18 x 24.5	Draeger®, Caleo® Incubators
1030091B	14 x 24	Ohio® IC (Round top)
TRANSPORT MATTRESSES		
1030094	13 x 26	Air Shields Isolette Models C86, C100, C200, CQT
1030098INT	10.75 x 23.4	Airborne® Transport Intubation Incubators
1030098	11.5 x 23.75	Airborne Transport Incubators Voyager and 185+ Models

All copyrights and registered trademarks are the property of their respective owners.



natus
pediatrics



Natus Medical Incorporated

1501 Industrial Road • San Carlos, CA 94070 USA | 1.800.303.0306 (US/CA) • Worldwide: 1.650.802.0400
www.natus.com