




- Sterile
- Single Use Needle Retracts and Locks After Use
- Two Controlled Depths of Penetration:
 - 1.0mm – yellow
 - 1.8mm – gray
- Latex-Free

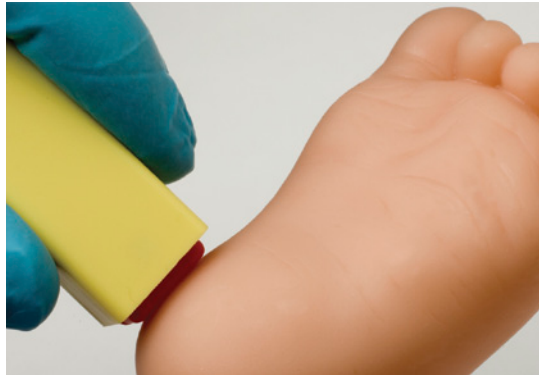

YELLOW FOR MICROPREEMIES TO FULL-TERM BABIES

GRAY FOR FULL-TERM BABIES AND ADULTS

Introducing SugarPlum – A glucose lancet for babies!

Wish you could draw just a single drop of blood for glucose tests? Now you can – with the SugarPlum glucose lancet.

- The same lancet style used for adult glucose tests – but sized for babies.
- Two controlled depths – 1.0mm depth for micropreemies to full-term babies; 1.8mm depth for full-term babies and adults.
- Economical.

SugarPlum & NeatNick™ – Perfect Companions For All Your Babies' Blood Draws!

Use SugarPlum for glucose tests that require just a drop of blood. Use our NeatNick sweeping blade lancet for PKU's, CBC's, and other tests. With NeatNick and SugarPlum you can choose the right lancet for all your babies' blood draws.

SugarPlum's patented design makes glucose draws quicker, easier and more comfortable for your babies. Press against the heel – SugarPlum automatically fires at the correct pressure. Before and after firing, the needle is safely concealed in the plastic body. 1.0mm depth for micropreemies to full-term babies (yellow); 1.8mm depth for full-term babies and adults (gray).

Ordering Information



For a product trial, questions, or the sales representative in your area, please call us.

ITEM #	DESCRIPTION	QUANTITY
<i>SugarPlum Glucose Lancet</i> Sterile, single-use disposable sharp. Latex-free.		
<i>For Micropreemies to Full-Term Babies</i>		
1040043	Un-bagged, 1.0mm depth (yellow)	100 per box
1040044	Individually Bagged, 1.0mm depth (yellow)	100 per box
<i>For Full-Term Babies and Adults</i>		
1040041	Un-bagged, 1.8mm depth (gray)	100 per box
1040042	Individually Bagged, 1.8mm depth (gray)	100 per box

PRODUCT SPECIFICATIONS:

1.0mm depth (g26) for Micropreemies to Full-Term Babies (yellow)
1.8mm depth (g23) for Full-Term Babies and Adults (gray)
SugarPlum is a sterile, single-use disposable "sharp" used to puncture an infant's heel for collecting a blood sample. Latex-free. See instructions for use.

Hawaii Medical
SugarPlum™

